

PowerSpin™ DXB Centrifuge

UNICO® is proud to present the PowerSpin™ DXB Centrifuge


With speeds up to 5500 RPM and horizontal spin rotor this centrifuge is a great instrument for a busy laboratory.

Featuring variable speed range of 500 – 5,500 RPM, to insure a clean and complete separation. Quiet and dependable, with a maintenance free motor, stainless steel tube shields it can be used to process blood, urine, fecal and other specimens.

The all digital front panel with 3 pre-programmed options

- 3 minutes at 4,000G (for BD Barricor™ tubes)
- 5 minutes at 3,200G
- 10 minutes at 2,000G

Custom settings are easy and intuitive with all digital panel and oversized display.

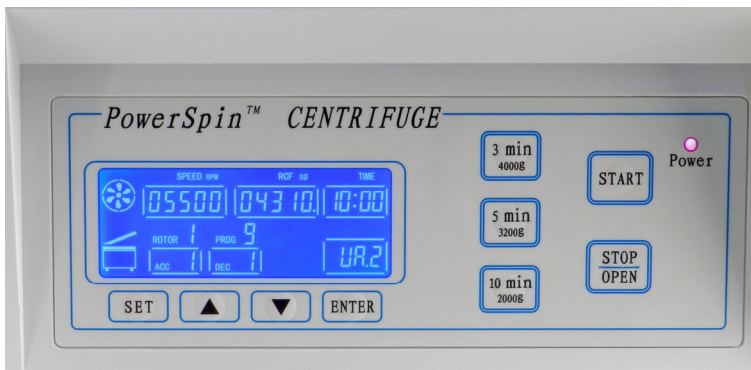

4000G for 3min is perfect for Barricor™ Tubes. BD Vacutainer® Barricor™ Plasma Blood Collection Tube are single-use, plastic, evacuated tubes with lithium heparin and a mechanical separator.

© 2016 BD. All rights reserved. Unless otherwise noted, BD, the BD Logo and all other trademarks are property of Becton, Dickinson and Company.

- Zero RPM Safety Latching System prevents lid from being opened while rotor is spinning
- Quiet running, only 50 decibels at 3500 RPM
- 99 minute digital timer, 1 sec. increment
- 500-5,500 RPM, 100 RPM increment

- 3 pre-programmed options for fast separation
- Remembers last setting for quick repeat runs
- Adjustable Acceleration and Deceleration
- Stainless steel rotor chamber for easy cleaning
- Available in 110V and 220V - Contact UNICO for details

MFG Part #	Description	UOM
C8763	PowerSpin DXB, 5500RPM Horizontal Spin Centrifuge. 110V	EA
C8763E	PowerSpin DXB, 5500 RPM Horizontal Spin Centrifuge. 220V	EA