


# Tube Adapter

Item: C800-34

**UNICO®**

**Directions for Use:** Hold a centrifuged Micro-Hematocrit Tube vertically against the chart with the bottom of the blood column at the 0% line. Slide the tube along the chart until the meniscus of the plasma crosses the 100% line. Read the height of the packed red cell column directly as percent cell volume.

## Micro-Hematocrit Capillary Tube Reader

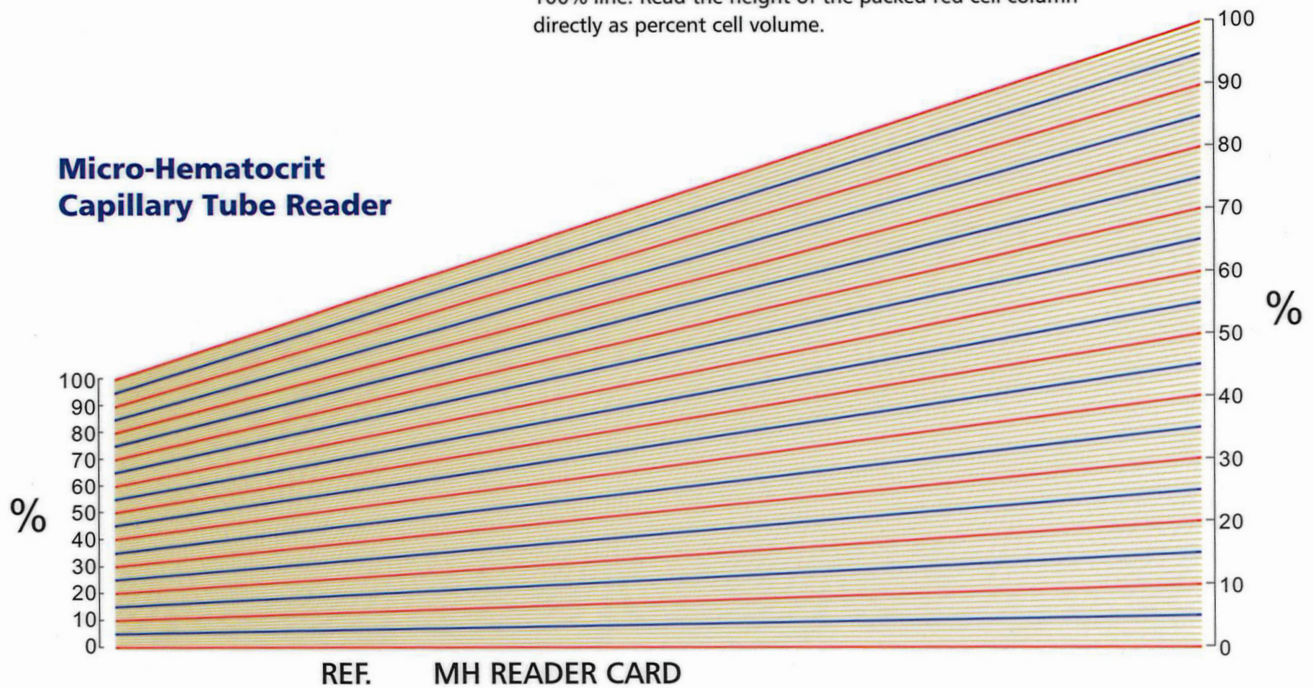

The C800-34 Microhematocrit Tube Adapter Set is specifically designed for use with the DX C8706 Centrifuge, offering a reliable solution for hematocrit testing. Each set includes two 4-place MH-tube adapters, which are capable of holding up to four microhematocrit tubes, ensuring proper balance and positioning during centrifugation. Additionally, the set comes with one Microhematocrit Reader Card, allowing for accurate and easy measurement of hematocrit values from the tubes after centrifugation. This adapter set is an essential tool for laboratories and clinical settings, providing efficiency and precision in hematocrit testing.

### COMPATIBILITY

**Centrifuge Model**

C8706

**Rotor Type**

6 Place Horizontal Rotor