


# PowerSpin™ DX Centrifuge

Model: C8704




- POWERSPIN™ DX CENTRIFUGE
- VARIABLE SPEED 1000-4000 RPM
- 4 PLACE HORIZONTAL ROTOR FOR 50 ML TUBES
- 1:00-99:99 MIN DIGITAL TIMER
- CAPACITY 4X50 ML TUBES
- ADAPTER PACKAGE AVAILABLE

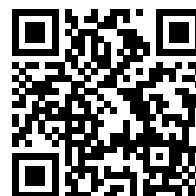

Boost your lab's efficiency with the C8704 DX Centrifuge. This robust benchtop centrifuge offers a 4-place swing-out rotor, providing ample capacity for 50ml tubes. Designed for convenience and reliability, the C8704 DX operates on standard 110V power, making it a practical addition to any laboratory setting. Experience quiet operation (50 decibels at 4000 RPM), precise control with the 99-minute digital timer and adjustable speed (500-4000 RPM), and time-saving features like last setting memory. The stainless steel rotor chamber ensures easy cleaning, and suction cup rubber feet minimize vibration. Safety is paramount with the Zero RPM Safety Latching System. For added versatility, you can order this unit with the adapter package which includes C8604-02 - Black Tube Adapter (set of 2), C8000-03-01 - Orange Short Tube Adapter (2), and C800-34 - Microhematocrit Tube Adapter with 2X4-place MH-Tube adapter, 2 tube shields, AND Reader Card. This adapter package expands the capabilities of your C8704 DX Centrifuge, allowing for a wider range of sample processing and enhancing your lab's overall productivity.



# C8704 Adapter Package

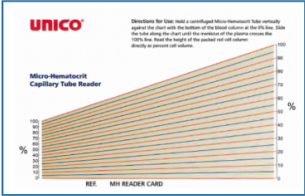

C8604-02 - Black Tube Adapter (set of 2)  
C8000-03-01 - Orange Short Tube Adapter (2)  
C800-34 - Microhematocrit Tube Adapter Set  
Microhematocrit Reader Card and Tools Included.



ADDITIONAL INFORMATION

TOP SELECTION	No	
ACCESSORIES	<b>Accessories (Sold Separately)</b>  C8600-01                      6 Place Horizontal Rotor, 2-15ml for DX Centrifuge C8600-02                      4 Place Horizontal rotor, 50ml DX Centrifuge  C8600-03M                      24 Place Fixed Angle Rotor, 2-15ml for DX Centrifuge, Metal (without Tube Shields)	
SPECIFICATION	<b>Speed (RPM)</b>	500-4000rpm
	<b>Maximum RCF (g)</b>	1848g
	<b>Rotor</b>	4 x 50ml Horizontal
	<b>Tube Size</b>	4 x 50ml
	<b>Motor</b>	DC
	<b>Power Input</b>	110V
SIZE AND WEIGHT	<b>Size</b>	13.75 inch Height (22 inch height with lid open) x 15 inch Deep x 11 inch Width
	<b>Shipping Size</b>	19"X18"X15"
	<b>Weight</b>	44 LBS
	<b>Shipping Weight</b>	44 LBS
WARRANTY	1 years	

Additional Options

Adapter Set	C8704 Adapter Package
-------------	-----------------------

