


# PowerSpin™ DX Centrifuge

Model: C8704E



- POWERSPIN™ DX CENTRIFUGE
- VARIABLE SPEED 1000-3400 RPM
- 4 PLACE HORIZONTAL ROTOR FOR 50 ML TUBES
- 1:00-99:99 MIN DIGITAL TIMER
- CAPACITY 4X50 ML TUBES
- POWER INPUT 220V

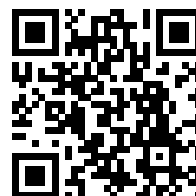

UNICO is proud to present the PowerSpin™ DX Centrifuge

With 3 rotor options this centrifuge is a great answer for every day laboratory tasks. Choose from 24 place fixed rotor for 2-15ml tubes (#C8724), 6 place horizontal rotor for 2-15ml tubes(#C8706) or 4 place horizontal rotor for up to 50ml tubes (#C8704)

The UNICO® PowerSpin™ DX C8700 Series centrifuge is engineered to meet the demanding needs of clinical, industrial, and biological labs. It offers reliable and efficient performance for blood separation and other material processing tasks. This centrifuge features digital variable speed control with a speed range from 500 RPM to 4000 RPM, an electronic timer with a range of 1:00 to 99:59, and a zero RPM safety lid latch that ensures safe operation by preventing the lid from opening while the unit is in motion. The intuitive controls allow users to select from preset programs for blood, urine, and fecal samples, or customize protocols easily using the set, up, down, and enter buttons. The unit's memory function retains the last time and speed settings, ensuring seamless operation for repetitive tasks. With a choice of rotors for standard 3 to 50 ml tubes, and optional adapters for smaller tubes, the C8700 Series is versatile enough for a wide range of laboratory applications. The centrifuge's maintenance-free brushless DC motor delivers quiet, reliable performance, while its rubber suction cup feet ensure stability during operation. Built with safety in mind, the unit will not run unless the lid is securely closed, and it features an automatic time and speed memory reset for convenience. Whether you're working in a clinical or research environment, the UNICO® PowerSpin™ DX C8700 Series offers a reliable and efficient solution for your centrifugation needs.





ADDITIONAL INFORMATION

TOP SELECTION	No	
ACCESSORIES	Accessories	
	C8600-01	6 Place Horizontal Rotor, 2-15ml for DX Centrifuge
	C8600-02	4 Place Horizontal rotor, 50ml DX Centrifuge
	C8600-03M	24 Place Fixed Angle Rotor, 2-15ml for DX Centrifuge, Metal (without Tube Shields)
	C800-02	Regular tube shield, pack of 2
	C800-03	Adapter for short tubes, pack of 10
	C800-04	Pediatric or microtube adapter, pack of 2
	C800-05	Tube Lift Adapters, Pk/12
SPECIFICATION	Speed (RPM)	500~4000 rpm
	Maximum RCF (g)	2508
	Rotor	4 place horizontal
	Tube Size	2-50ml *
	Motor	DC
	Power Input	220V
SIZE AND WEIGHT	Size	13.75 inch Height (22 inch height with lid open) x 15 inch Deep x 11 inch Width
	Shipping Size	19"X18"X15"
	Weight	31 LBS
	Shipping Weight	36 LBS
WARRANTY	1 years	

