


PowerSpin™ DX Centrifuge

Model: C8763H

- POWERSPIN™ DXB CENTRIFUGE
- VARIABLE SPEED 500-5500 RPM
- 3 PRE-PROGRAMMED OPTIONS FOR BD BARRICOR™ AND OTHER TUBES
- 6 PLACE HORIZONTAL ROTOR
- 60 MIN DIGITAL TIMER
- CAPACITY 6X15 ML TUBES
- 24 PLACE TUBE HOLDSTER™ RACK
- POWER INPUT 110V

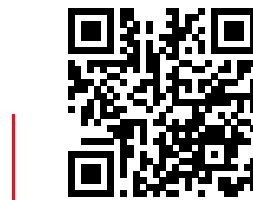

UNICO® is proud to present the PowerSpin™ DXB Centrifuge. With speeds up to 5500 RPM and horizontal spin rotor this centrifuge is a great instrument for a busy laboratory

Featuring variable speed range of 500 – 5,500 RPM, to insure a clean and complete separation. Quiet and dependable, with a maintenance free motor, stainless steel tube shields it can be used to process blood, urine, fecal and other specimens.

The all digital front panel with 3 pre-programmed options

3 minutes at 4,000G (for BD Barricor™ tube)

1800 RPM for 5 minutes for Urine

10 minutes at 2,000G

Custom settings are easy and intuitive with all digital panel and oversized display.

4000G for 3min is perfect for Barricor™ Tubes. BD Vacutainer® Barricor™ Plasma Blood Collection Tube are single-use, plastic, evacuated tubes with lithium heparin and a mechanical separator.

© 2016 BD. All rights reserved. Unless otherwise noted, BD, the BD Logo and all other trademarks are property of



ADDITIONAL INFORMATION

ACCESSORIES	Accessories <ul style="list-style-type: none"> C8600-01 6 Place Horizontal Rotor, 2-15ml for DX Centrifuge C8600-02 4 Place Horizontal rotor, 50ml DX Centrifuge C8600-03M 24 Place Fixed Angle Rotor, 2-15ml for DX Centrifuge, Metal (without Tube Shields) C800-02 Regular tube shield, pack of 2 C800-03 Adapter for short tubes, pack of 10 C800-04 Pediatric or microtube adapter, pack of 2 C800-05 Tube Lift Adapters, Pk/12
SPECIFICATION	<ul style="list-style-type: none"> Speed (RPM) 500~4000 rpm Maximum RCF (g) 4000 Rotor 6 place horizontal Tube Size 2-15ml * Motor DC Power Input 220V
SIZE AND WEIGHT	<ul style="list-style-type: none"> Size 13.75 inch Height (22 inch height with lid open) x 15 inch Deep x 11 inch Width Shipping Size 19"X18"X15" Weight 31 LBS Shipping Weight 36 LBS

