


# PowerSpin™ DX Centrifuge

Model: C8764H

- POWERSPIN™ DXB CENTRIFUGE
- VARIABLE SPEED 500-5500 RPM
- 3 PRE-PROGRAMMED OPTIONS FOR BD BARRICOR™ AND OTHER TUBES
- 4 PLACE HORIZONTAL ROTOR FOR 50 ML TUBES
- 60 MIN DIGITAL TIMER
- CAPACITY 4X50 ML TUBES
- 24 PLACE TUBE HOLDSTER™ RACK
- POWER INPUT 110V

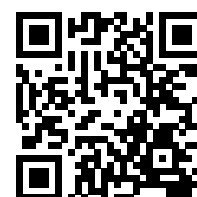

UNICO® is proud to present the PowerSpin™ DXB Centrifuge. With speeds up to 5500 RPM and horizontal spin rotor this centrifuge is a great instrument for a busy laboratory

Featuring variable speed range of 500 – 5,500 RPM, to insure a clean and complete separation. Quiet and dependable, with a maintenance free motor, stainless steel tube shields it can be used to process blood, urine, fecal and other specimens.

The all digital front panel with 3 pre-programmed options

3 minutes at 4,000G (for BD Barricor™ tube)

1800 RPM for 5 minutes for Urine

10 minutes at 2,000G

Custom settings are easy and intuitive with all digital panel and oversized display.

4000G for 3min is perfect for Barricor™ Tubes. BD Vacutainer® Barricor™ Plasma Blood Collection Tube are single-use, plastic, evacuated tubes with lithium heparin and a mechanical separator.

© 2016 BD. All rights reserved. Unless otherwise noted, BD, the BD Logo and all other trademarks are property of



**ADDITIONAL INFORMATION**

<p><b>ACCESSORIES</b></p>	<p><b>Accessories</b></p> <p>C8600-01                      6 Place Horizontal Rotor, 2-15ml for DX Centrifuge</p> <p>C8600-02                      4 Place Horizontal rotor, 50ml DX Centrifuge</p> <p>C8600-03M                    24 Place Fixed Angle Rotor, 2-15ml for DX Centrifuge, Metal (without Tube Shields)</p> <p>C800-02                        Regular tube shield, pack of 2</p> <p>C800-03                        Adapter for short tubes, pack of 10</p> <p>C800-04                        Pediatric or microtube adapter, pack of 2</p> <p>C800-05                        Tube Lift Adapters, Pk/12</p>
<p><b>SIZE AND WEIGHT</b></p>	<p><b>Size</b>                            13.75 inch Height (22 inch height with lid open) x 15 inch Deep x 11 inch Width</p> <p><b>Shipping Size</b>                19"X18"X15"</p> <p><b>Weight</b>                         31 LBS</p> <p><b>Shipping Weight</b>             36 LBS</p>

