

PowerSpin™ BX Centrifuge

Model: C882E



- POWERSPIN™ BX CENTRIFUGE
- VARIABLE SPEED 1000-11000 RPM
- 24 PLACE MICROHEMATOCRIT ROTOR
- 99 MIN DIGITAL TIMER
- CAPACITY 24 MICROHEMATOCRIT TUBES
- POWER INPUT 220 V




UNICO is proud to introduce the latest addition to the PowerSpin™ Family – the BX model high speed centrifuge with 3 rotor options:

8-place Blood Tube rotor, 24 place Microhematocrit and 24 place Micro-Tube rotor

The state of the art design incorporates advanced safety features like Zero RPM Lid Lock and Auto Unbalance Detection, Rotor Recognition and Shutdown.

The built in oversized LCD Display provides real time data to the user. With the push of a button you can speed up or slow down the rotor within 1,000 - 11,000 RPM range (Max 5000 RPM with the Blood Tube Rotor). Maintenance free brushless DC motor is powerful and super quiet. The cover of the Microhematocrit rotor incorporates a reader, allowing for instant reading and eliminating the need for additional specimen handling.

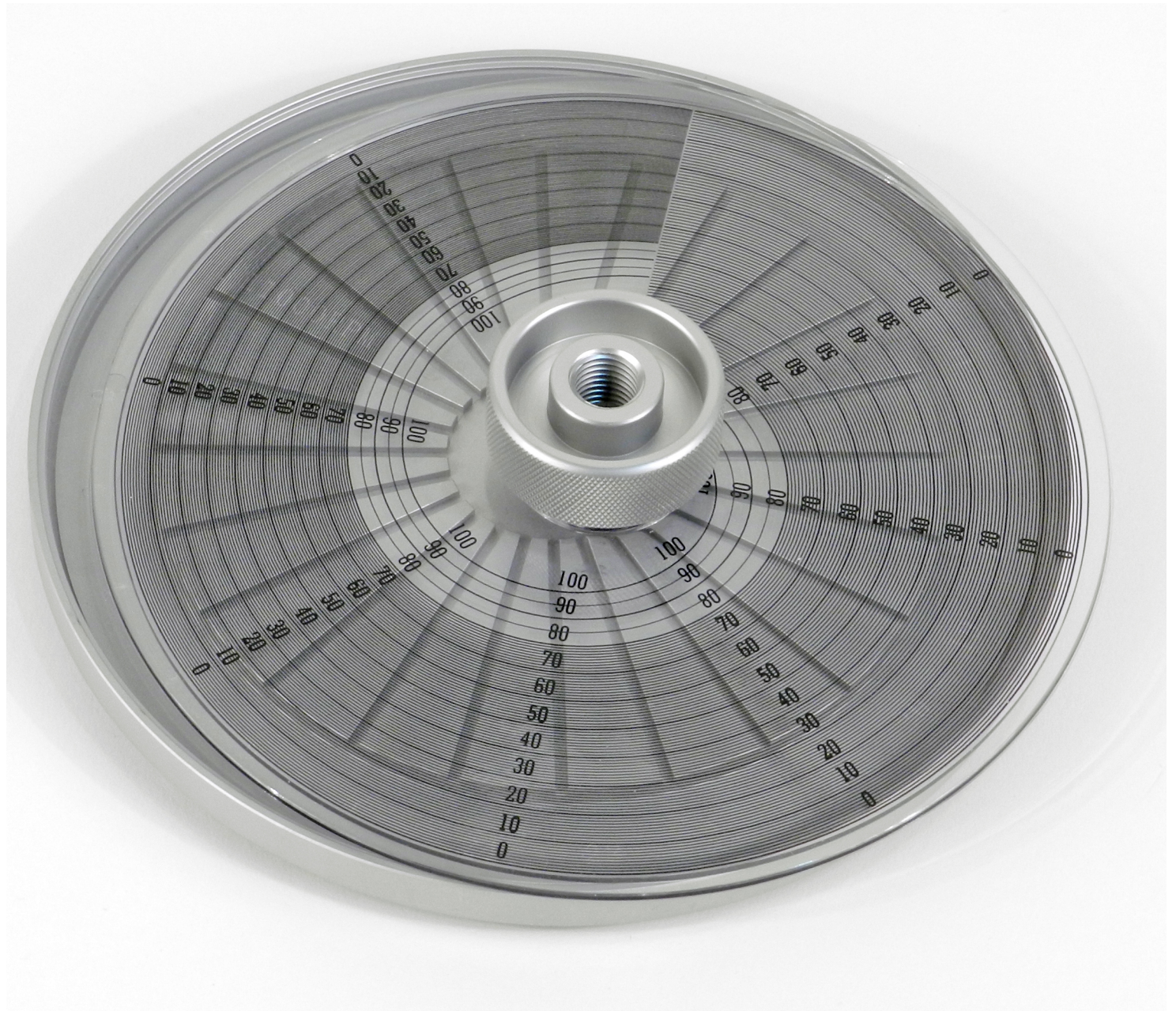

ADDITIONAL INFORMATION

TOP SELECTION	No	
ACCESSORIES	Accessories	
	C880-08	Fixed Angle Rotor, 8 place
	C880-02	Microhematocrit Rotor, 24 place
	C880-03	Microtube rotor, 24 place
	C880-10	BX Rotor Rack
	C881-03	Short tube adapter, for 13x75mm tubes
	C881-04	Micro / Pediatric tube adapter for 0.5-2.0ml tubes
SPECIFICATION	Rotor	Fixed Angle, 24 place
	Speed Range	1000 - 11000
	Maximum RCF(g)	11634
	Nominal Capacity	1.5 x 75 mm
	Maximum Test Tube Size	1.5 x 75 mm
	Power Input	220V
SIZE AND WEIGHT	Size	10.5 in x 12.5 in x 9 in (270 mm x 320 mm x 220 mm)
	Shipping Size	14 in x 15 in x 12 in (350 mm x 375 mm x 300 mm)
	Weight	15.5 lbs (C881 & C882); 16.5 lbs (C883)
	Shipping Weight	19 lbs (C881 and C882) ; 20 lbs (C883)
WARRANTY	1 years	

