

PowerSpin™ BX Centrifuge


Model: C884



- POWERSPIN™ BX CENTRIFUGE
- VARIABLE SPEED 1000-11000 RPM
- 8 PLACE FIXED ANGLE AND 24PLACE MICROHEMATOCRIT ROTOR
- 99 MIN DIGITAL TIMER
- POWER INPUT 110V



The BX Centrifuge is a versatile benchtop centrifuge designed for a range of laboratory applications. This particular model offers exceptional flexibility with two included rotors: a larger 8-place rotor designed for 10mL tubes (8 x 10mL capacity) ideal for larger volume sample processing, and a 24-place microhematocrit rotor for efficient and accurate hematocrit determinations. This dual-rotor capability eliminates the need for purchasing separate centrifuges, saving valuable lab space and budget. The BX Centrifuge is powered by a standard 110V power supply, making it readily adaptable to most laboratory environments. The built in oversized LCD Display provides real time data to the user. With the push of a button you can speed up or slow down the rotor within 1,000 - 11,000 RPM range (Max 5000 RPM with the Blood Tube Rotor). Maintenance free brushless DC motor is powerful and super quiet. The cover of the Microhematocrit rotor incorporates a reader, allowing for instant reading and eliminating the need for additional specimen handling.

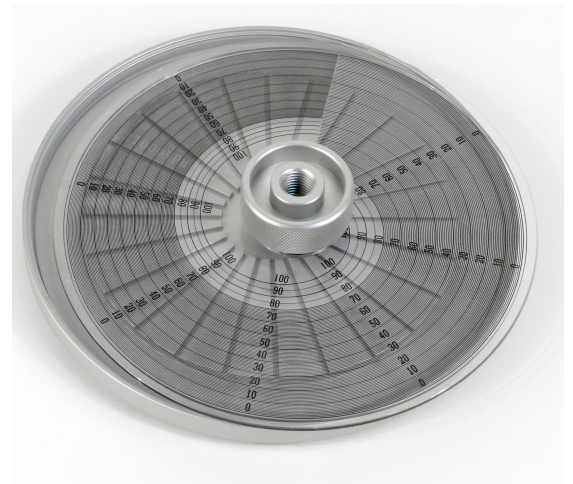

ADDITIONAL INFORMATION

TOP SELECTION	No	
ACCESSORIES	Accessories (Sold Separately) C880-08 Fixed Angle Rotor, 8 place C880-02 Microhematocrit Rotor, 24 place C880-03 Microtube rotor, 24 place C880-10 BX Rotor Rack C881-03 Short tube adapter, for 13x75mm tubes C881-04 Micro / Pediatric tube adapter for 0.5-2.0ml tubes	
SPECIFICATION	Rotor	Microcentrifuge, 24 place
	Speed Range	1000 - 11000
	Maximum RCF(g)	11093
	Nominal Capacity	24 x 0.5-2.0 ml
	Maximum Test Tube Size	5 x 10mm
SIZE AND WEIGHT	Size	10.5 in x 12.5 in x 9 in (270 mm x 320 mm x 220 mm)
	Shipping Size	14 in x 15 in x 12 in (350 mm x 375 mm x 300 mm)
	Weight	15.5 lbs (C881 & C882); 16.5 lbs (C883)
	Shipping Weight	19 lbs (C881 and C882) ; 20 lbs (C883)
WARRANTY	1 years	

