


PowerSpin™ BX Centrifuge

Model: C884E

- POWERSPIN™ BX CENTRIFUGE
- VARIABLE SPEED 1000-11000 RPM
- 8 PLACE FIXED ANGLE AND 24PLACE MICROHEMATOCRIT ROTOR
- 99 MIN DIGITAL TIMER
- POWER INPUT 220V

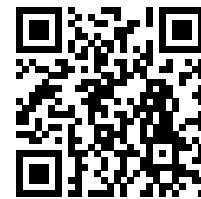


UNICO is proud to introduce the latest addition to the PowerSpin™ Family – the BX model high speed centrifuge with 3 rotor options:

8-place Blood Tube rotor, 24 place Microhematocrit and 24 place Micro-Tube rotor

The state of the art design incorporates advanced safety features like Zero RPM Lid Lock and Auto Unbalance Detection, Rotor Recognition and Shutdown.

The built in oversized LCD Display provides real time data to the user. With the push of a button you can speed up or slow down the rotor within 1,000 - 11,000 RPM range (Max 5000 RPM with the Blood Tube Rotor). Maintenance free brushless DC motor is powerful and super quiet. The cover of the Microhematocrit rotor incorporates a reader, allowing for instant reading and eliminating the need for additional specimen handling.

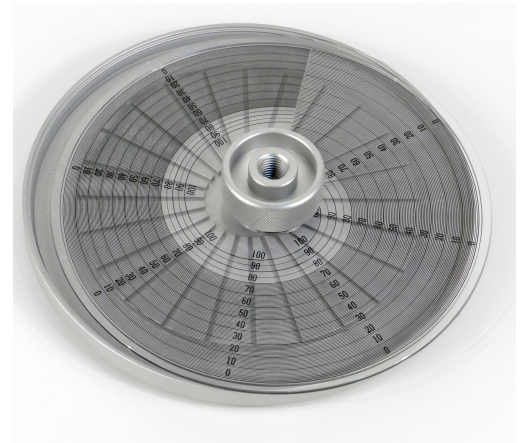

ADDITIONAL INFORMATION

ACCESSORIES	<p>Accessories</p> <p>C880-08 Fixed Angle Rotor, 8 place</p> <p>C880-02 Microhematocrit Rotor, 24 place</p> <p>C880-03 Microtube rotor, 24 place</p> <p>C880-10 BX Rotor Rack</p> <p>C881-03 Short tube adapter, for 13x75mm tubes</p> <p>C881-04 Micro / Pediatric tube adapter for 0.5-2.0ml tubes</p>
SPECIFICATION	<p>Rotor Microcentrifuge, 24 place</p> <p>Speed Range 1000 - 11000</p> <p>Maximum RCF(g) 11093</p> <p>Nominal Capacity 24 x 0.5-2.0 ml</p> <p>Maximum Test Tube Size 5 x 10mm</p> <p>Power Input 220V</p>
SIZE AND WEIGHT	<p>Size 10.5 in x 12.5 in x 9 in (270 mm x 320 mm x 220 mm)</p> <p>Shipping Size 14 in x 15 in x 12 in (350 mm x 375 mm x 300 mm)</p> <p>Weight 15.5 lbs (C881 & C882); 16.5 lbs (C883)</p> <p>Shipping Weight 19 lbs (C881 and C882) ; 20 lbs (C883)</p>

