


PowerSpin™ MH Centrifuge

Model: CMH30E




- POWERSPIN™ MH MICRO-HEMATOCRIT CENTRIFUGE
- FIXED SPEED 12000 RPM
- 24 PLACE MICROHEMATOCRIT ROTOR
- 5 MIN MECHANICAL TIMER
- POWER INPUT 220V
- MICROHEMATOCRIT READER INCLUDED.

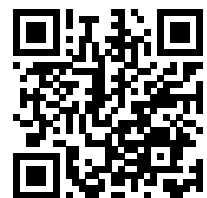

Model C-MH30 is made of ABS plastic with inner metal chamber. A safety switch built into the handle assembly disconnects the power to the motor whenever the latch is lifted.

The 24 place microhematocrit rotor is an aluminum shell for durability and lighter weight

The brushed motor is mounted on rubber mounts to provide quieter running and less vibration. cleaned and disinfected.

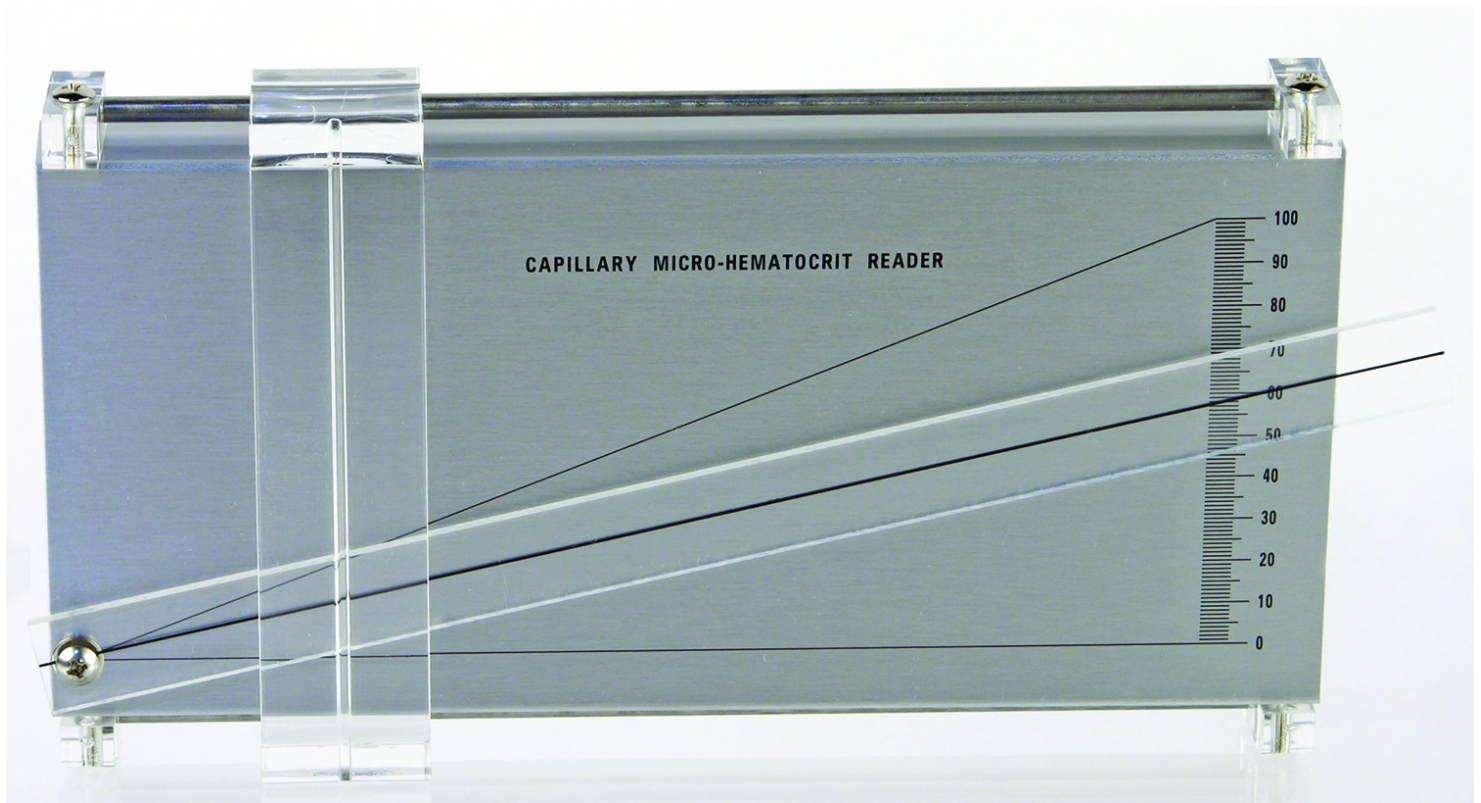

ADDITIONAL INFORMATION

TOP SELECTION	No	
ACCESSORIES	Accessories	
	CMH30-Reader	Microhematocrit reader
	CMH30-03	Gasket Replacement For CMH30 Centrifuge
SPECIFICATION	Rotor	24 place
	Power Input	220V
SIZE AND WEIGHT	Size	11.6 in x 10.3 in x 9.3 in (295 x 261 x 235 mm)
	Shipping Size	12.7 in x 12.7 in x 12.2 in (323 x 323 x 311 mm)
	Weight	17.4 lbs (7.9 kg)
	Shipping Weight	19.4 lbs (8.8 kg)
WARRANTY	1 years	

