


# Tube Rocker

Model: L-TTR200DE

- UNICO TEST TUBE ROCKER
- UP TO 30 TUBES/VIALS CAPACITY
- EASY OPERATION TOP PAUSE BUTTON
- FIXED SPEED OF 25 ROCKING CYCLES PER MINUTE
- POWER INPUT 220 V

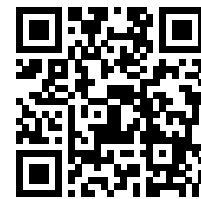

UNICO Rock-IT's now feature a PAUSE button, a useful feature designed to allow the user the opportunity to pause the rocking motion while adding or removing tubes. The feature also extends the life of the rocker motor by eliminating the need to hold the platform to add/remove tubes.

UNICO Rock-IT Tube Rockers are heavy-duty, yet nearly silent mixers designed for years of everyday use. Smooth and gentle movement at 25 RPM assures a uniform, thoroughly mixed specimen, while preserving the integrity of the fragile blood cells.

Get your hematology controls ready to use in no time for accurate results. The reversible non-slip pad is autoclavable and will accommodate larger and hard to mix tubes. Rock-IT's are an economical choice for your lab; let the Rock IT work for you!



**ADDITIONAL INFORMATION**

<b>SPECIFICATION</b>	<b>Capacity</b>	30 tubes
	<b>Speed</b>	25cpm
	<b>Voltage</b>	220V
<b>SIZE AND WEIGHT</b>	<b>Size</b>	8.6 in x 4.7 in x 5.0 in
	<b>Shipping Size</b>	10.6 in x 6.7 in x 6.9 in
	<b>Weight</b>	4.0 lbs (1.8 kg)
	<b>Shipping Weight</b>	5.0 lbs (2.3 kg)

