


Tube Mixer

Model: MTR22

- UNICO MULTI-PURPOSE MIXER
- DIGITAL PROGRAMMABLE OPTIONS
- 360° ROTATION AND VORTEXING COMBINATION
- VARIABLE SPEED AND TIMER
- 2 PLATFORMS FOR VARIOUS SIZE TUBES AND VIALS
- 90 - 240 VOLT UNIVERSAL POWER SUPPLY

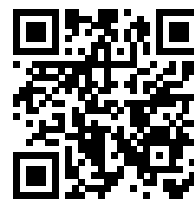

UNICO MTR22 Multi-Mix has the capability to perform three modes of mixing, all of which are programmable, and can all be used in the same "run". Or, if you just want to rotationally mix, set the unit to rotate. It is very simple to operate and program, runs can be programmed in seconds. You can mix many different sizes of tubes ranging from 1.5 ml to 50 ml with the two interchangeable mixing platforms that come standard with the unit. Tubes are held in-place with unique rubber flexible grippers; that surround and firmly hold tubes of varying diameters.



ADDITIONAL INFORMATION

SIZE AND WEIGHT	Size	13.5 in x 7 in x 7.5 in (342 x177 x 192 mm)
	Shipping Size	15 in x 8 in x 8 in (375 x 285 x 285 mm)
	Weight	6.7 lbs (3 kg)
	Shipping Weight	7.5 lbs (3.8 kg)

