


Clinical Laboratory Incubator

Model: L-CU200E

- UNICO CLINICAL LABORATORY INCUBATOR
- 20L CHAMBER CAPACITY
- 2 REMOVABLE SHALVES
- DOUBLE DOOR DESIGN
- TEMPERATURE AMBIENT TO 70° C
- POWER INPUT 220 VOLT

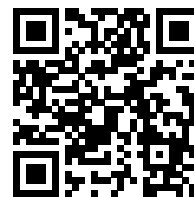

UNICO incubators are designed for plates, tubes (slants), urine and general microbiology applications. The spacious stainless steel-lined interiors can also be used as fluid warmers. These incubators feature a double-door design, solid exterior door and a see-thru inner acrylic door.



ADDITIONAL INFORMATION

SPECIFICATION	Chamber Capacity	20 Liters (0.7 cubic feet)
	Chamber Size	Internal Dimensions 12.3"W x 12.1"D x 9"H
	Calibration	Digital control & thermometer access port
	Temperature Range	Ambient to 70°C (Ambient to 155° F)
	Shelves	Two removable shelves, perforated for air flow
	Temperature Control	Digital Control
	Display	Digital Display
	Power/Consumption	1.4 Amps / 150 Watts
	Heating	Convection
	Stability	Will hold constant temperature +/- 0.5° C
	Power input	Available in 110V and 220V

