


# Infinity Phase Microscope

Model: IP755

- SEIDENTOPF BINOCULAR HEAD
- WF10X PAIRED EYEPIECES
- 10X,40XR,100XR DIN INFINITY PLAN PHASE OBJECTIVES
- BUILT-IN MECHANICAL STAGE WITH PHASE CONDENSER
- HALOGEN ILLUMINATION

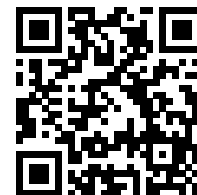

The UNICO IP755 Infinity Plan Phase Binocular Microscope is a high-performance tool designed for advanced microscopy applications. This versatile microscope is ideal for research, clinical, and educational environments where detailed phase contrast imaging and exceptional optical clarity are essential.

Featuring WF10X widefield eyepieces, the IP755 provides a broad and comfortable field of view, enhancing user experience and reducing eye strain during extended use. The microscope is equipped with a set of infinity plan objectives, including a 4X brightfield objective and 10X, 40X, and 100X (oil immersion) phase contrast objectives. These objectives are engineered to deliver precise, distortion-free images and excellent phase contrast, making it easier to observe and analyze transparent specimens.

The IP755 includes a 6V20W halogen illuminator with adjustable intensity control, allowing users to tailor the lighting conditions for optimal contrast and visibility. This feature ensures consistent and customizable illumination, improving the overall effectiveness of the phase contrast technique.

With its advanced optical system, versatile phase contrast objectives, and user-friendly design, the UNICO IP755 Infinity Plan Phase Binocular Microscope offers outstanding performance and reliability. It is an excellent choice for professionals seeking high-quality imaging and precise specimen analysis in a range of microscopy applications.



ADDITIONAL INFORMATION

TOP SELECTION	No
ACCESSORIES	<b>Bodies</b>
	SEIDENTOPF BINOCULAR BODY WITHOUT EYEPIECES FOR IP730 SERIES MICROSCOPE
	IP730-3004
	SEIDENTOPF TRINOCULAR BODY, 30 DEGREE INCLINED WITHOUT EYEPIECES FOR IP730 SERIES MICROSCOPE
	IP730-3005
	SEIDENTOPF BINOCULAR BODY WITHOUT EYEPIECES FOR IP50 SERIES MICROSCOPE
	IP750-3004
	SEIDENTOPF TRINOCULAR BODY, 30 DEGREE INCLINED WITHOUT EYEPIECES FOR IP750 SERIES MICROSCOPE
	IP750-3005
	<b>Optics</b>
	10X EXTRA WIDEFIELD EYEPIECE FN 22MM FOR IP730 AND IP750 MICROSCOPE
	IP750-1502
	4X DIN ACHROMAT OBJECTIVE, N.A. 0.10
	IP730-2101
	10X DIN ACHROMAT OBJECTIVE, N.A. 0.25
	IP730-2102
	40X DIN ACHROMAT OBJECTIVE, N.A. 0.65, RETRACTABLE FRONT LENS
	IP730-2103
	100X DIN ACHROMAT OBJECTIVE, N.A. 1.25, RETRACTABLE FRONT LENS
	IP730-2104
	PHASE SET W/ PLAN PHASE 10X-40X-100X OBJECTIVES, TURRET CONDENSER WITH "0" POSITION
	IP730-8030
	4X INFINITY ACHROMATIC OBJECTIVE FOR IP750 INFINITY MICROSCOPE
	IP750-2101
	10X INFINITY ACHROMATIC OBJECTIVE FOR IP750 INFINITY MICROSCOPE
	IP750-2102
	40X INFINITY ACHROMATIC OBJECTIVE FOR IP750 INFINITY MICROSCOPE
	IP750-2103
	100X INFINITY ACHROMATIC OBJECTIVE FOR IP750 INFINITY MICROSCOPE
	IP750-2104
	4X INFINITY PLAN ACHROMATIC OBJECTIVE FOR IP750 INFINITY MICROSCOPE
	IP750-2301



	IP750-2302	10X INFINITY PLAN ACHROMATIC OBJECTIVE FOR IP750 INFINITY MICROSCOPE
	IP750-2303	40X INFINITY PLAN ACHROMATIC OBJECTIVE FOR IP750 INFINITY MICROSCOPE
	IP750-2304	100X INFINITY PLAN ACHROMATIC OBJECTIVE FOR IP750 INFINITY MICROSCOPE
	IP750-2422	INFINITY PLAN PHASE OBJECTIVE 10X FOR IP750 SERIES
	IP750-2423	INFINITY PLAN PHASE OBJECTIVE 40X FOR IP750 SERIES
	IP750-2424	INFINITY PLAN PHASE OBJECTIVE 100X FOR IP750 SERIES
	IP750-8030	PHASE SET W/ INFINITY PLAN PHASE 10X-40X-100X OBJECTIVES, TURRET CONDENSER WITH "0" POSITION
	<b>Miscellaneous Accessories</b>	
	IP750-5001	MECHANICAL STAGE WITHOUT ABBE CONDENSER FOR IP730 AND IP750 MICROSCOPE
	IP750-6002	PHASE CONTRAST CONDENSER, TURRET, NA1.25 WITH 10X-40X-100X PHASE RING AND "0" BRIGHT FIELD POSITION
<b>SPECIFICATION</b>	IP750-6003	DARKFIELD CONDERNSER FOR IP730 AND IP750 SERIES MICROSCOPE
	IP730-8002	DUST COVER FOR IP730 AND 750 SERIES MICROSCOPE
	IP750-7006	6V 20W HALOGEN BULB, 4MM PIN FOR IP730 AND IP750 MICROSCOPE
	IP750-8901	HARD CARRYING CASE FOR IP730 AND IP750 SERIES MICROSCOPE
<b>SIZE AND WEIGHT</b>	<b>Head</b>	Binocular
	<b>Eyepiece</b>	WF10X
	<b>Objectives</b>	4x Plan 10x 40XR 100XR Infinity Phase Plan
	<b>Stage</b>	MECH.
	<b>Condenser</b>	Abbe N.A.1.25
	<b>Illuminator</b>	6V/20W /Halogen
	<b>Size</b>	10 in x 8.5 in x 15.4 in
	<b>Shipping Size</b>	11 in x 14 in x 18 in
	<b>Weight</b>	18 lbs





	<b>Shipping Weight</b>	20 lbs
--	------------------------	--------

